

Art Plus Family

Premium technology for **superb quality**



Premium technology for superb quality

FONA Art Plus is an advanced digital panoramic system with the latest premium CdTe-CMOS sensor technology. It combines extraordinary image quality, sharpness and contrast with advanced diagnostic tools. Due to the special sensor technology and continuous autofocus function, FONA Art Plus always delivers properly focused images, in both Pan and Pan Ceph versions.

CdTe-CMOS sensor

CCD/CMOS sensor

X-ravs




CdTe-CMOS Sensor (X-ray directly converted to electrical signal)



Premium Art Plus technology

Direct conversion of X-rays to electrical signals increases the image clarity and contrast while reducing the readout time and dose.

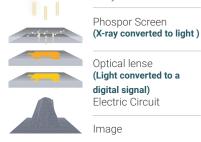


Traditional technology

Traditional acquisition uses conversion of X-rays to light before detection, causing some information loss, blurring and readout delays.

Wide focal zone up to 30 mm


Due to patented algorithms, the focal one is up to 30 mm wide, instead of traditional 4 mm. Your region of interest will always be in focus.

Latest sensor and autofocus technology for the best image quality

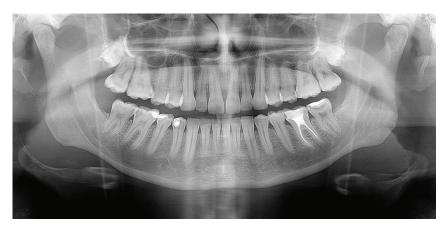





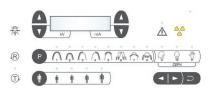

Only 2 steps to go

Simply select program (8 Pan or 3 Ceph) and patient size (4 options). All exposure settings will be set automatically. Presets can also be adjusted manually.




Single sensor for both applications Get the best quality panoramic and cephalometric images. Just click the sensor to the position where you need it.



•	1 3	
Patient size p	oresets	

(1) Large. (2) Medium. (3) Small. (4) Child.

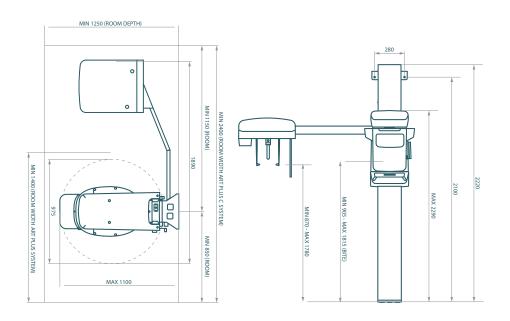

Italian Radiology

Through more than 20 years of experience in production, we design and develop extraoral and intraoral equipments for dental clinics, radiographic centers, hospitals and research institutes. Advanced technology connected with your devices for analysis and diagnosis. Based in Milano, Italy, we have more than 10.000 panoramic units and 60.000 intraorals installed worldwide. We provide with international support and assistance always caring about our customers and people.

Latest technology and people caring

Digital dentistry

Intraoral imaging

Extraoral imaging

Intraoral sensors



Based in Italy, available worldwide



Copyright FONA s.r.l., All rights reserved, Subject to technical changes and errors in text, BRO002, v2, ENG, 202301

Count on us









HotLine Support



FONA s.r.l. Via G.Galilei 11 20057 Assago (MI) - Italy

Contacts

info@fonaitaly.com www.fonadental.com





Clinical Cases

Online Training

Remote Support